

C8-cAMP caproKit™ 50 Reactions

List of components

C8-cAMP Capture Compound™

The C8-cAMP caproKit™ allows a selective isolation of cAMP-binding proteins. The synthetic C8-cAMP Capture Compound™ (Figure 1) uses cyclic adenosyl-monophosphate (cAMP) as selectivity function to interrogate native proteins. cAMP is attached via an aminohexylamino group at C8 to the Capture Compound Scaffold. The Capture Compound Mass Spectrometry (CCMS) technology enables analysis, discovery and characterization of cAMP-binding proteins through an efficient reduction of the complexity of the proteome.

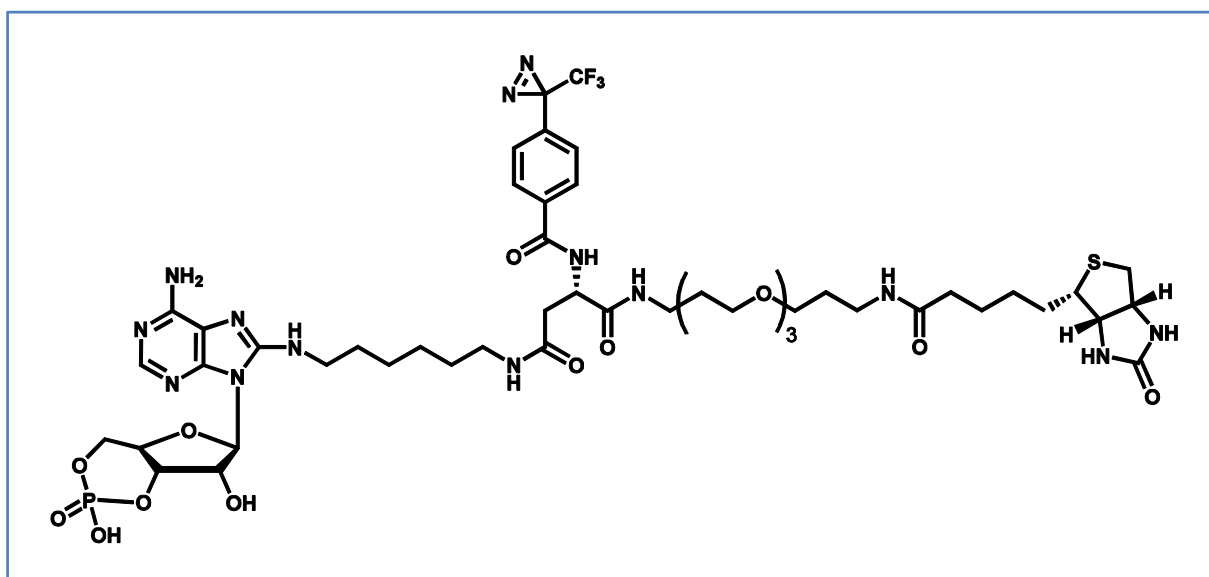


Figure 1: C8-cAMP Capture Compound™ for selective isolation of cAMP-binding proteins using cAMP as selectivity function.

Item No	Component	Amount	Buffer composition
2-3010-050	Streptavidin coated magnetic beads (10 mg/ml SA-MB)	2.8 ml	Dynal Dynabeads MyOne™ Streptavidin C1 (Invitrogen)
2-2100-050	Capture buffer 1 (5x CB1)	1.2 ml	HEPES, KOAc, Mg(OAc) ₂ , Glycerol, Triton X-100, pH 7.5
2-2200-050	Wash buffer 1 (5x WB1)	25 ml	Tris-HCl, EDTA, NaCl, Octyl-β-D-glucopyranoside, pH 7.9
2-1030-050	C8-cAMP Capture Compound™ (C8-cAMP-CC, 100 μM)	1.4 ml	Water
2-4030-050	C8-cAMP competitor (4 mM)	1.5 ml	HEPES, pH 7.5
2-5040-050	PKA RI (43 kDa, 16 μM)	5 μl	NaCl, K ₂ HPO ₄ , EDTA, Glycerol, β-mercaptoethanol, pH 6.8
3-4011-000	12 PCR Tube strips 0.2 ml (AB-1114)	5	

Note: The PKA RI solution will not freeze due to the glycerol in the storage buffer. It is recommended to aliquot the C8-cAMP Capture Compound™ solution upon arrival to avoid multiple freezing/thawing cycles of the stock solution. Protect the C8-cAMP Capture Compound™ from direct light.


Separate 1.5 ml 5x WB1 in a fresh 2.0 ml tube for the assay. Dilute the rest of the 5x WB1 in a ratio of 1:5 in aqua bidest for washing steps and store at -20 °C to -18 °C. **Do not freeze the Streptavidin coated magnetic beads!** All solutions must be entirely thawed and mixed before usage.

Storage notification

Item No	Component	During shipment (max 3 days)	After receipt
2-3010-050	Streptavidin coated magnetic beads (10 mg/ml SA-MB)	4 to 8 °C	4 to 8 °C
2-2100-050	Capture buffer 1 (5x CB1)	4 to 8 °C	- 20 to -18 °C
2-2200-050	Wash buffer 1 (5x WB1)	4 to 8 °C	- 20 to -18 °C
2-1030-050	C8-cAMP Capture Compound™ (C8-cAMP-CC, 100 μM)	4 to 8 °C	- 20 to -18 °C
2-4030-050	C8-cAMP competitor (4 mM)	4 to 8 °C	- 20 to -18 °C
2-5040-050	PKA RI (43 kDa, 16 μM)	4 to 8 °C	- 20 to -18 °C

Specified Functionality

Significant band (SDS-PAGE/silver stain) with 0.5 μg (11 pmol) of PKA RI and significant competition with C8-cAMP, when the capturing protocol described in the C8-cAMP caproKit™ guideline is applied, caproBox™ and only kit components are used.

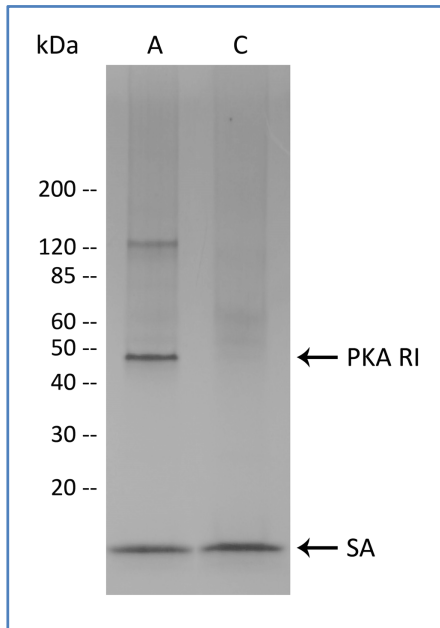

A: Capture assay with PKA RI

C: Control of "A" using 1 mM C8-cAMP as competitor

SA: Streptavidin from Streptavidin coated magnetic beads


Figure 2: Capture assay (A) and C8-cAMP competition control (C) of the positive control enzyme PKA RI analyzed by SDS-PAGE/silver stain.

Stability

The C8-cAMP caproKit™ is stable under storage conditions for 6 months. After first use microbial contamination may occur.

Please read the material safety data sheet for this product at www.caprotec.com

Berlin, 08/23/2010



Head of Quality Control

Contact and order information:

Headquarters

caprotec bioanalytics GmbH

Volmerstrasse 5
D-12489 Berlin

Phone: +49 30 63 92 39 90

Fax: +49 30 63 92 39 89

Web: www.caprotec.com

Email: sales@caprotec.com

caprotec Inc., USA

15 New England Executive Office Park
Burlington, MA 01803, USA

Phone: +1 781 685 4992

Fax: +1 781 685 4601

Web: www.caprotec.com

Email: sales@caprotec.com

Ordering information: caprotec part number: 1-1030-050

Copyright

© 2008-2010 caprotec bioanalytics GmbH. All rights reserved. Reproduction in whole or in part only with permission of caprotec bioanalytics GmbH.

Trademarks

caproKit, caproBox, caprotec, caproBeads, caproMag, ImproMed and Capture Compound are trademarks of caprotec bioanalytics GmbH. All other used tradenames or trademarks belong to their respective owners.

Products & Services

CCMS technology is made available as ready to use caproKit reagents and services.

For more information please visit www.caprotec.com

Or contact us. Email: info@caprotec.com

Phone: +49 30 6392 4004

Products and Services are for Research use only.

© 2008-2010 caprotec bioanalytics GmbH